


KIT COMPONENTS

Component	Amount
96-well plate precoated with anti-human TNF β antibody	1 Plate
Protein Standard: Lyophilized recombinant human TNF β	2 tubes, 10 ng/tube
Sample Diluent Buffer	30 ml
Biotinylated Antibody (Anti-human TNF β)	130 μ l (100x)
Antibody Diluent Buffer	12ml
Avidin-Biotin-Peroxidase Complex (ABC) Solution	130 μ l (100x)
ABC Diluent Buffer	12 ml
Tetramethylbenzidine (TMB) Color Developing Agent	10 ml
TMB Stop Solution	10 ml

Storage

Store at 4°C. Cell Applications, Inc. recommends using the kit within 6 months of order.

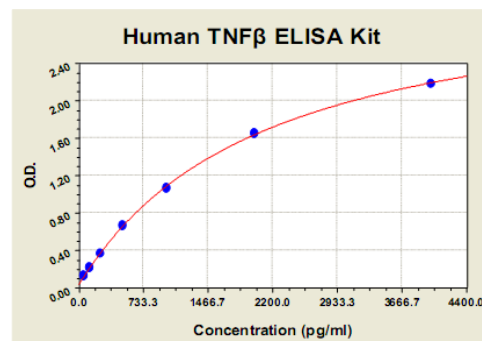

X	pg/ml	0.0	62.5	125	250	500	1000	2000	4000
Y	O.D.	0.029	0.134	0.225	0.369	0.668	1.064	1.656	2.180

Figure 1: TNF β Standard Curve. Using the Human TNF β ELISA Kit, O.D. data was graphed against TNF β protein concentration. The TMB reaction was incubated at 37°C for 15 min.

BACKGROUND

The human TNF-beta, also known as lymphotoxin-A (LTA), is an important cytokine involved in the development of secondary lymphoid organs and inflammatory responses. This factor is produced predominantly by mitogen-stimulated T-lymphocytes and leukocytes. The factor is secreted also by fibroblasts, astrocytes, myeloma cells, endothelial cells, epithelial cells and a number of transformed cell lines. The synthesis of TNF-beta is stimulated by interferons and IL2. Some pre-B-cell lines and Abelson murine leukemia virus-transformed pre-B-cell lines constitutively produce TNF-beta. The human TNF-beta protein contains 205 amino acids. There are two forms of TNF-beta: cell membrane and extracellular soluble form which is in blood stream and biological fluids. The soluble form of TNF-beta is usually a homotrimer with a relative molecular mass of 60 to 70 kDa, whereas the membrane form of TNF-beta is a heteromeric complex with lymphotoxin beta (TNF β , LT β , TNFSF3). The human TNF-beta shares 35% identity and 50% homology in amino acid sequence with the human TNF-alpha. The biological function of TNF-beta is mediated largely by TNF-alpha receptor 1 and TNF-alpha receptor 2. Several studies suggested that TNF-beta can also recognize LIGHT (TNFSF14) receptor. The polymorphism of TNF-beta gene in either the coding region or the promoter region has been associated with Crohn disease and myocardial infarction.¹

TNF-beta acts on a plethora of different cells. In general, TNF-beta and TNF-alpha display similar spectra of biological activities *in vitro* systems, although TNF-beta is often less potent or displays apparent partial agonist activity. TNF-beta is cytolytic or cytostatic for many tumor cells. Hemorrhagic necrosis of tumors induced by TNF-beta *in vivo* is probably the result of an inhibition of the growth of endothelial cells and the activity of TNF-beta as an anti-angiogenesis factor.² In monocytes TNF-beta induces the terminal differentiation and the synthesis of G-CSF. TNF-beta is a mitogen for B-lymphocytes. In neutrophils TNF-beta induces the production of reactive oxygen species. It is also a chemoattractant for these cells, increases phagocytosis, and also increases adhesion to the endothelium. In addition, TNF-beta induces the synthesis of GM-CSF, G-CSF, IL1, collagenase, and prostaglandin E2 in fibroblasts. TNF-beta promotes the proliferation of fibroblasts and is involved probably in processes of wound healing *in vivo*. TNF-beta inhibits the growth of osteoclasts and keratinocytes.³

References

- Ruddle, N.H.: Curr. Opin. Immunol. 4: 327-32, 1992
- Aggarwal, B.B. & Singh, S.: Growth Fact. Cytok. Health Dis. 2: 727-75, 1997
- Aggarwal, B.B.: Immunol Ser. 56:61-78, 1992

ELISA OVERVIEW

Cell Applications ELISA Kits are based on standard sandwich enzyme-linked immunosorbent assay technology. Freshly prepared standards, samples, and solutions are recommended for best results.

- Prepare test samples.
- Prepare a protein standard of the target protein.
- Add test samples and standards to the pre-coated 96-well plate. Do not wash.
- Add biotinylated detection antibodies. Wash.
- Add Avidin-Biotin-Peroxidase Complex (ABC) Solution. Wash.
- Add Tetramethylbenzidine (TMB) Color Developing Agent, containing HRP substrate.
- Add TMB Stop Solution
- Subject the plate to analysis.

NOTES:

- Before using the kit, quick spin tubes to bring down all solution to the bottom of tube.
- Duplicate assay wells are recommended for both standard and sample testing.
- Do not let the 96-well plate dry, this will lead to inactivation of plate components.
- When diluting samples and reagents, ensure that they are mixed completely and evenly.
- Pre-warm diluted ABC and TMB solutions at 37°C for 30 min before use to avoid variable temperature effects.
- For washes, use TBS or PBS. Do not touch well walls.
- A protein standard is included in the kit. A protein standard detection curve should be generated with each experiment, no more than 2 hours prior to the experiment.
- The user will determine sample dilution fold by estimation of target protein amount in samples.

FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR CLINICAL USE.



ELISA PROTOCOL

Preparation of Test Samples


- Process Test Samples in the following manner:
 - Cell culture supernate, tissue lysate or body fluids:** Remove particulates by centrifugation.
 - Serum:** Allow the serum to clot in a serum separator tube (about 30 min) at room temperature. Centrifuge at approximately 1000 x g for 15 min. Analyze the serum immediately or aliquot and store frozen at -20° C.
 - Plasma:** Collect plasma using heparin, EDTA as an anticoagulant. Centrifuge for 15 min at 1000 x g within 30 min of collection. Analyze immediately or aliquot and store frozen at -20° C. Citrate is not recommended as the anticoagulant.
- Estimate the concentration of the target protein in the sample and select a proper dilution factor such that the diluted target protein concentration falls within the 62.5-4000 pg/ml standard curve range. Depending on the sample, several trial dilutions may be necessary. Dilute the sample using the provided diluent buffer, mixing well. Suggested working dilutions of samples are as follows:

Target Protein Concentration Range	Sample Working Dilution	Sample Vol.	Diluent Buffer Vol.
40-400 ng/ml	1:100	1 μ l	99 μ l
4-40 ng/ml	1:10	10 μ l	90 μ l
62.5-4000 pg/ml	1:2	50 μ l	50 μ l
\leq 62.5 pg/ml	n/a	100 μ l	n/a

- If samples will be assayed within 24 hours, store at 2-8°C. For long-term storage, aliquot and freeze samples at -20°C. Avoid repeated freeze-thaw cycles.

Preparation of Standard Solutions (62.5-4000 pg/ml)

- Reconstitute the Lypophilized Recombinant Protein to make a 10,000 pg/ml TNF β solution. Add 1 ml Sample Diluent Buffer to a tube of lypophilized protein, keep the tube at room temperature for 10 min. Mix thoroughly.
- Aliquot 0.6 ml of the Sample Diluent Buffer to an eppendorf tube, and label as 4000 pg/ml Protein Standard.
- Add 0.4 ml of the mixed 10,000 pg/ml TNF β solution to the eppendorf tube containing 0.6 ml diluent buffer and mix to make a 4000 pg/ml TNF β solution.
- Label 6 eppendorf tubes with the TNF β protein concentrations to be prepared by serial dilution: 2000pg/ml, 1000pg/ml, 500pg/ml, 250pg/ml, 125pg/ml, 62.5pg/ml.
- Aliquot 0.3 ml of the Sample Diluent Buffer to the labeled tubes.
- Serially dilute the protein standards into their respectively labeled tubes. Transfer 0.3 ml from the 4000pg/ml TNF β Solution to the 2000pg/ml eppendorf tube and mix thoroughly. Transfer 0.3 ml of the 2000 pg/ml solution to the 1000pg/ml tube and mix thoroughly, Transfer 0.3 ml of the 1000pg/ml solution to the 500pg/ml tube and mix, and so on to make the 250, 125, and 62.5 pg/ml solutions.
- Store at 4°C until use.

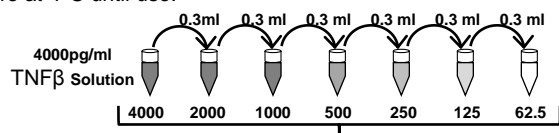

Loading the 96-well Plate TNF β Protein Standards [pg/ml]

- Aliquot 0.1 ml of the sample diluent buffer into a control well to serve as the Blank. This will yield the O.D.450(Blank) reading.
- Aliquot 0.1 ml of the standard solutions of the **Preparation of Standard Solutions** (62.5-4000 pg/ml) into empty wells of the precoated 96-well plate. Duplicate measurements of standards are recommended.

- Aliquot 0.1 ml of each properly diluted test sample to empty wells prepared in **Step 2**. Duplicate measurements of each test sample are recommended.
- Cover the 96-well plate and incubate at 37°C for 90 min.
- During the **Step 14** incubation period, prepare a stock of Biotinylated 1:100 Antibody Working Solution. Count the number of reactions and multiply by 0.1 ml/well for the Working Solution total volume (preparation of 1-2 reactions in excess of the number of wells is recommended). Dilute the Biotinylated Antibody to 1:100 in Antibody Diluent Buffer and mix thoroughly. Use the working solution within 2 hours.
- Upon completion of the 90 min incubation of **Step 14**, remove the cover of the 96 well plate and discard plate well contents. Blot the plate onto paper towels or other absorbent material. DO NOT let the wells completely dry at any time.
- Add 0.1 ml of the Biotinylated 1:100 Antibody Working Solution (prepared in **Step 15**) to each well and incubate the plate at 37°C for 60 min.
- During the incubation period of **Step 17**, prepare a stock of ABC Working Solution. Count the number of reactions and multiply by 0.1 ml/well for the Working Solution total volume (preparation of 1-2 reactions in excess of the number of wells is recommended). Dilute the ABC Stock Solution to 1:100 in ABC Diluent Buffer and mix thoroughly. Pre-warm the ABC working solution at 37°C for 30 min before use. Use the working solution within 1 hour.
- Upon completion of the 60 min incubation of **Step 17**, wash the plate 3 times with 0.3 ml TBS or PBS. For each wash, leave washing buffer in the wells for 1-2 min. Discard the washing buffer and blot the plate onto paper towels or other absorbent material.
- Add 0.1 ml of prepared ABC Working Solution (prepared in **Step 18**) to each well and incubate the plate at 37°C for 30 min.
- During the incubation period of **Step 20**, pre-warm TMB Color Developing Agent at 37°C for 30 min before use.
- Upon completion of the 30 min incubation of **Step 20**, wash the plate 5 times with 0.3 ml TBS or PBS. For each wash, leave the washing buffer in the wells for 1-2 min. Discard the washing buffer and blot the plate onto paper towels or other absorbent material.
- Add 90 μ l of the pre-warmed TMB Color Developing Agent into each well and incubate at 37°C for 10-15 min (shades of blue can be seen in the wells with the four most concentrated Protein Standard Solutions; the other control wells should show no obvious color).
- Add 0.1 ml of the TMB Stop Solution to each well. The acidic stop solution will change the mixture color to yellow. The yellow intensity is proportional to the amount of target protein captured by the plate.
- Read the O.D. absorbance at 450nm in a microplate reader within 30 min after adding the stop solution. These readings are the O.D.450(Reading).

Calculating Protein Concentration

- For all wells, determine O.D.450(Relative):

$$\text{O.D.450(Relative)} = \text{O.D.450(Reading)} - \text{O.D.450(Blank)}$$
- Plot the standard curve: plot O.D.450(Relative) of each standard solution (Y) vs. the respective concentration of the standard solution (X). See **Figure 1** for a typical standard curve.
- The TNF β concentration of the samples can be interpolated from the standard curve. Multiply the interpolated concentration by the dilution factor to obtain the target protein concentration in the sample.

FOR RESEARCH USE ONLY. NOT FOR DIAGNOSTIC OR CLINICAL USE.

